



Medik is a group focusing on surgical room solution and nursing care solution
Medik is a group focusing on surgical room solution and nursing care solution. We started to manufacture surgical table & lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.

Qualified with ISO9001, ISO13485, CE, Medik has successfully distributed the surgical products covering surgical table & lamp & pendant to international hospitals both domestic and overseas including many countries in Europe.

info@Mediksurgical.com | +86 139 6223 1292

MK-D620620W

LED Surgical Light




MK-D620620W

LED SURGICAL LIGHT

FEATURES

- Multi-Segment Lamp Arm, provides enhanced flexibility and stability
- High-Precision Sensor guarantees stable lighting in complex environments
- Endoscope Eye Protection Mode offers multiple options suitable for different surgical procedures
- Adjustable Color Temperature allows customization of the lighting color temperature to meet specific surgical requirements and preferences
- Intuitive interface of LED Control Panel for easy lighting adjustments, ensuring efficiency.
- Elegant Design enhances operating room aesthetics with a beautiful exterior design
- Shadow Compensation reduces glare for improved visibility during surgeries.



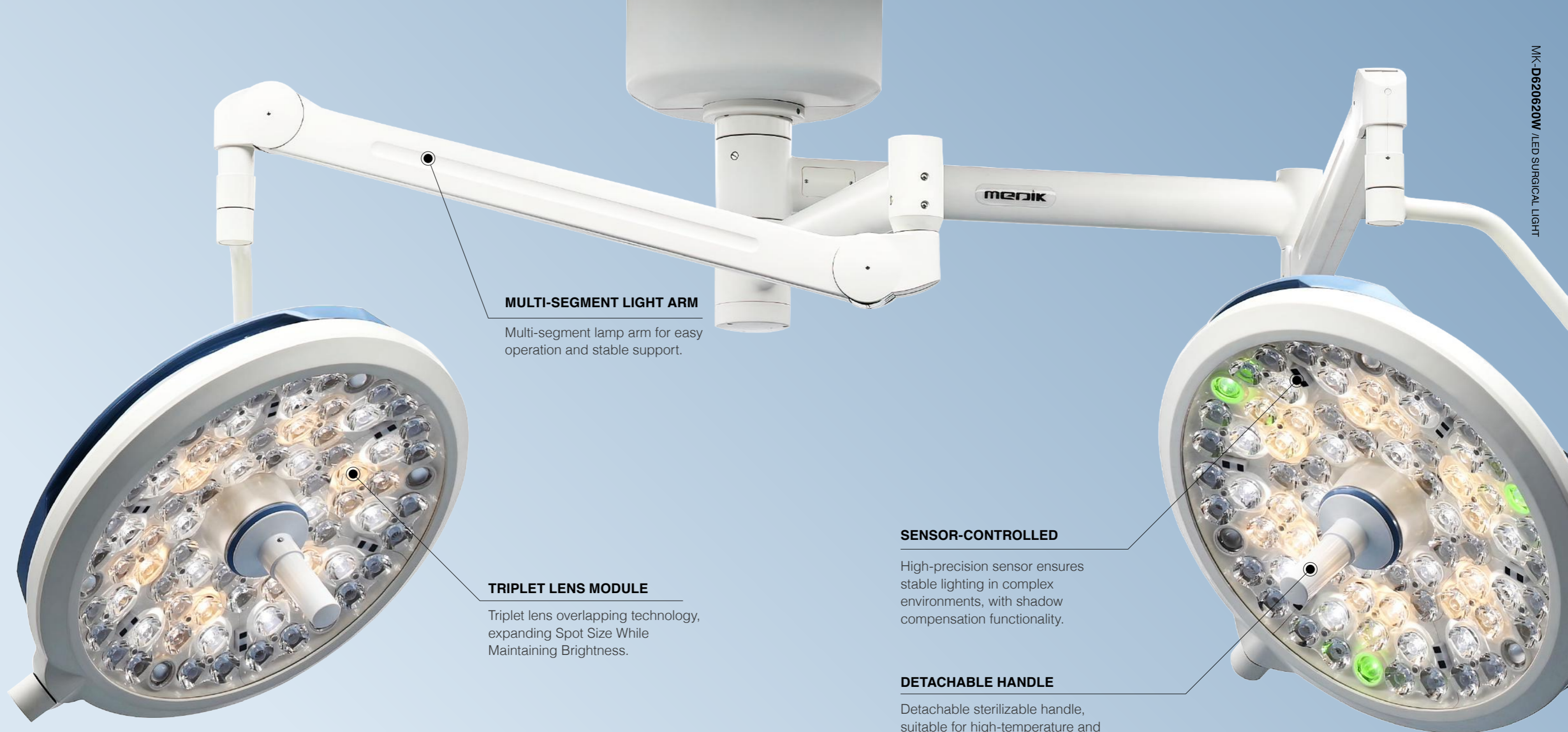

MULTI-SEGMENT LIGHT ARM

Multi-segment lamp arm for easy operation and stable support.

TRIPLET LENS MODULE

Triplet lens overlapping technology, expanding Spot Size While Maintaining Brightness.

SENSOR-CONTROLLED

High-precision sensor ensures stable lighting in complex environments, with shadow compensation functionality.

DETACHABLE HANDLE

Detachable sterilizable handle, suitable for high-temperature and high-pressure disinfection

MK-D620620W

LED SURGICAL LIGHT

The mk-d620620w led surgical light sets a new standard in operating room illumination. With its striking design and advanced features, it delivers unparalleled performance for surgical procedures. The led control panel allows for precise light adjustment, while shade compensation ensures optimal visibility. Designed to meet various surgical needs, this wall-mounted light combines beauty and functionality seamlessly.

MK-D620620W

LED SURGICAL LIGHT



Luminance(Lux)	40,000-160,000
Diameter of light head(mm)	620&620mm
Number of bulbs	78&78
Color temperature(K)	3000-5000 adjustable
Color render index(Ra)	≥98%
Spot Diameter (mm)	120-350
LED bulb brand	Germany Osram
LED Power (VA)	100
Illuminance depth(mm)	1500
Bulb service life(H)	>60000
Voltage	AC110-240V 50/60H

ERGONOMIC HANDLE

Streamlined handle design for enhanced aesthetics and ergonomic comfort.




THIN LIGHT HOLDER

Optimized with a thin light holder for fast cooling, energy efficiency, and extended lifespan.



ADJUSTABLE COLOR TEMPERATURE

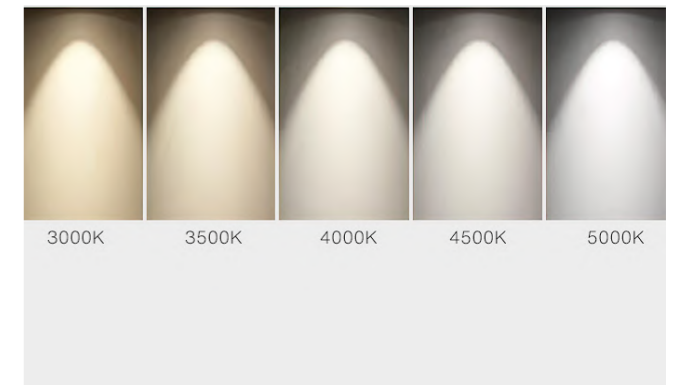

LED DISPLAY CONTROL

Convenient and Precise Brightness Adjustment.

