

MK-D700JL

Mobile Surgical Light (Flower Type)



*Better care
Better comfort!*

MK-D700JL

MOBILE SURGICAL LIGHT (FLOWER TYPE)

MK-D700JL Flower LED surgical light offers optimal lighting with four calibrated color temperatures. Users can choose warm orange tones or cooler warm white, adapting to visual preferences. The latest LED technology ensures a user-friendly configuration, making it a perfect lighting solution for various surgeries.

FEATURES

- Set individually from 3,800 K to 5,500 K for high-contrast light tailored to various tissue types.
- With Endoscopy Mode, Switch to White LEDs for optimal conditions during endoscopic interventions.
- LAS (Light Adjustment System) technology maintains 100% lighting power even after years of use.
- Facilitates easy disinfection, enhancing overall hygiene and infection control.
- Easily positioned slim and curved stand to get close to the lighting area.
- Superior optics with adjustable color temperature.
- Better shadow-free performance.
- Maneuverable control of light heads to reduce infection possibilities.
- Excellent heat management and wound care.
- Single, dual, ceiling, wall, and mobile mounting options

ITEM	PARAMETER
Model	MK-D700JL
Light Intensity At (1M)	60,000-160,000Lux
LED Bulb Quantity	61PCS
Bulb Life (HRS)	≧50,000
Color Temp (Kelvin)	3800 ~5500K (4Steps Adjustable)
LED Bulb Size	20MM
Color Rendering Index (Ra)	≧96
Color rendering index R9(Red)	98
Color Temperature Adjustable	3800-5500K
Total radiant flux density	364W/m ²
Rotating Angle Of Arm	>540°

PRODUCT DETAILS

LED Bulb

Powered by long-lasting LED bulbs, ensures an extended lifespan and minimal heat emission.




Removable adjustable handle

Removable handle for easy disinfection, essential for surgical lights in operating rooms.



Multi-segmented arms

Featuring multi-segmented arms allows flexible adjustment at various angles. Tailor the illumination precisely to your needs, offering adaptability for diverse applications.

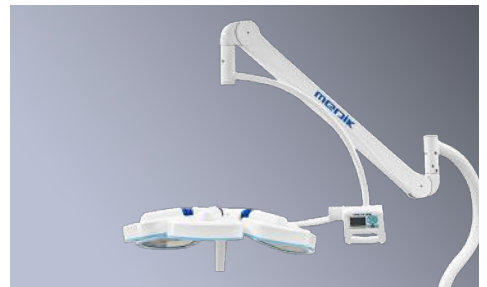

Clear switch

LCD screen for precise control in a compact design for enhanced user experience.



PRODUCT DETAILS

Mobile Base

With a mobile base that allows easy movement and positioning of the illumination light.



Luminaire housing

Flower LED ensures hygiene with a flow-optimized luminaire housing designed for laminar airflow ceilings, providing an IP54 protection level.



ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

www.mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US