

MK-E250ZL

LED Examination Light



*Better care
Better comfort!*

MK-E250ZL

LED EXAMINATION LIGHT

DESCRIPTION

The MK-E250ZL is a LED Examination Light equipped with German OSRAM light beads, It ensures optimal brightness and clarity for medical procedures and examinations,The slim profile of the light not only adds a sleek aesthetic to the examination room but also maximizes maneuverability.

MAIN FEATURES

1. The light-head case comes in two colors, it is unique.
2. Constant current drive design to extend its working life.
3. Its LED Bulb is Energy saving and has low power consumption.
4. German Osram LED bulb, it can work for more than 50000hours.
5. Overall aluminum with integrated plastic shell, better heat dissipation effect.
6. Maximize lens light effect, adjustable high and low color temperature, with illumination memory function.

TECHNICAL PARAMETER

Model	MK-E250ZL
Voltage	AC110-240V 50/60HZ
Spot diameter	180-220mm (adjustable)
Illuminance	>40000Lx (4 levels adjustable)
Color temperature	4000-5000k (adjustable)
Color rendering index	≥90% (adjustable)
LED lamp bead life	≥50000h

PRODUCT DETAILS

Handles

It incorporates a handle that is easy to adjust around the light allowing doctors to position the light precisely.



Ultra-thin light holder

The slim light head enhances the look and efficiently dissipates heat of light.



Clear adjustment icon

Featuring a clear adjustment icon, this examination light offers intuitive control for doctors.



Basement

The C-shaped base provides enhanced stability for added security.



PRODUCT DETAILS


OSRAM light beads

With the Osram light beads which known for exceptional brightness and clarity, ensure a long lifespan of over 50,000 hours, delivering enduring performance.



Adjustable head

The light head tilts up 90 degrees, down 90 degrees, and rotates horizontally up to 300 degrees for versatile positioning.

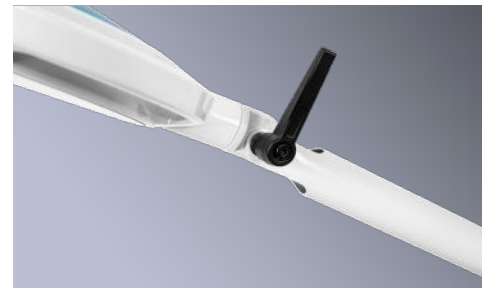

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales. The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

www.mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US