

YA-03E

Electric Surgical Table



*Better care
Better comfort!*

YA-03E

ELECTRIC SURGICAL TABLE

DESCRIPTION

YA-03E is an electric surgical table with comprehensive functionality and high performance. Hand remote controller with embedded railing controller provides double security during operation. T shape base provides medical staff a closer access to patient with high stability.

FEATURES

- Divided into 6 parts which can meet the needs of various surgeries
- All metal parts are made of high-grade 304 stainless steel, which is corrosion-resistant and easy to clean
- T-shaped base, which prevents tipping and allows doctors to get as close to the patient as possible, making surgery easier
- Equipped with remote control and column control system
- The transparent platform and sliding functions are suitable for C-arm
- Kidney bridge with adjustable height
- 4 wheels, which are very stable after locking

STANDARD SPECIFICATION

Table Length	2070mm (81.5")
Table Width	550mm (21.5")
Table Height Range	700~1000mm (27.5"~39.3")
Head section up/down	45° / 90°
Back section up/down	75° / 20°
Leg section up/down	15° / 90°
Leg section outward	90°
Trendelenburg	25°
Reverse trendelenburg	25°
Lateral Tilt (left/right)	15° / 15°
kidney bridge	120mm (4.7")
Sliding	300mm (11.8")
Loading capacity	250kg (550lb)

Lowest height position
700mm



Highest height position
1000mm



Back plate up
0-75°



Back plate down
0-20°



Reverse trendelenburg
0-25°



Trendelenburg
0-25°



Lateral tilt left
0-15°



Lateral tilt right
0-15°



Horizontal sliding
300mm



Head section up
0-45°



Head section down
0-90°



Leg section down
0-90°



Leg section outward
0-90°



Kidney elevator
120mm



Leg section up
0-15°





Kidney bridge

built-in kidney bridge, can be adjusted up to 120mm through the crank



Leg section adjustment device

Legs are adjustable via knobs and gas springs




Head section adjustment device

Head adjusted via mechanical gears



Armrest adjustment device

Angle of armrest adjusted via gas spring

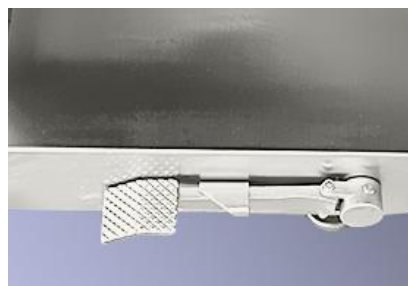

Castor

Moves via 4 castors, stable and does not shake after locking



Hand remote hanger



Column control system



Hand remote

ABOUT US

Medik is a group focusing on surgical room solution and nursing care solution

Medik is a group focusing on surgical room solution and nursing care solution. We started to manufacture surgical table & lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.

Qualified with ISO9001, ISO13485, CE, Medik has successfully distributed the surgical products covering surgical table & lamp & pendant to international hospitals both domestic and overseas including many countries in Europe.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US