

YA-B104

Orthopedic Traction Extension



*Better care
Better comfort!*

YA-B104

ORTHOPEDIC TRACTION EXTENSION

Made of stainless steel and with double joints, the YA-B104 orthopedic traction extension offers the best performance in terms of durability and adjustment range. It is suitable for C-arms and can provide doctors with the largest field of view due to the air suspending design. The multi-dimensional and multi-angle adjustments of various parts can assist doctors in obtaining the best surgical position.

MAIN FEATURES

- Main part is made of stainless steel, antibacterial, durable, anti-corrosion.
- Dual joint design ,larger adjustment range and wider surgical field of view.
- The boot can be adjusted in various heights, angles and distances to provide the best surgical position.
- The traction mechanism has a precise scale display.
- Equipped with an auxiliary trolley, it is easy to install and use, and is suitable for various operating tables.


TECHNICAL PARAMETER

| | |
|---|--|
| Height adjustment | Can follow the height of surgical table |
| Carbon Fiber Pole Outreach | 360° |
| Traction stroke | ≥200mm |
| Horizontal movement range of knee support | ≥550mm |
| Vertical angle adjustment range of traction mechanism | 50°~140° |
| Overall size | 1800*350mm (Width and height can be adjusted to the operating table) |

PRODUCT DETAILS

Dual Joint

Double-joint traction frame, with a wider adjustable range, enables more surgical positions and provides the doctor with the largest surgical field of view.

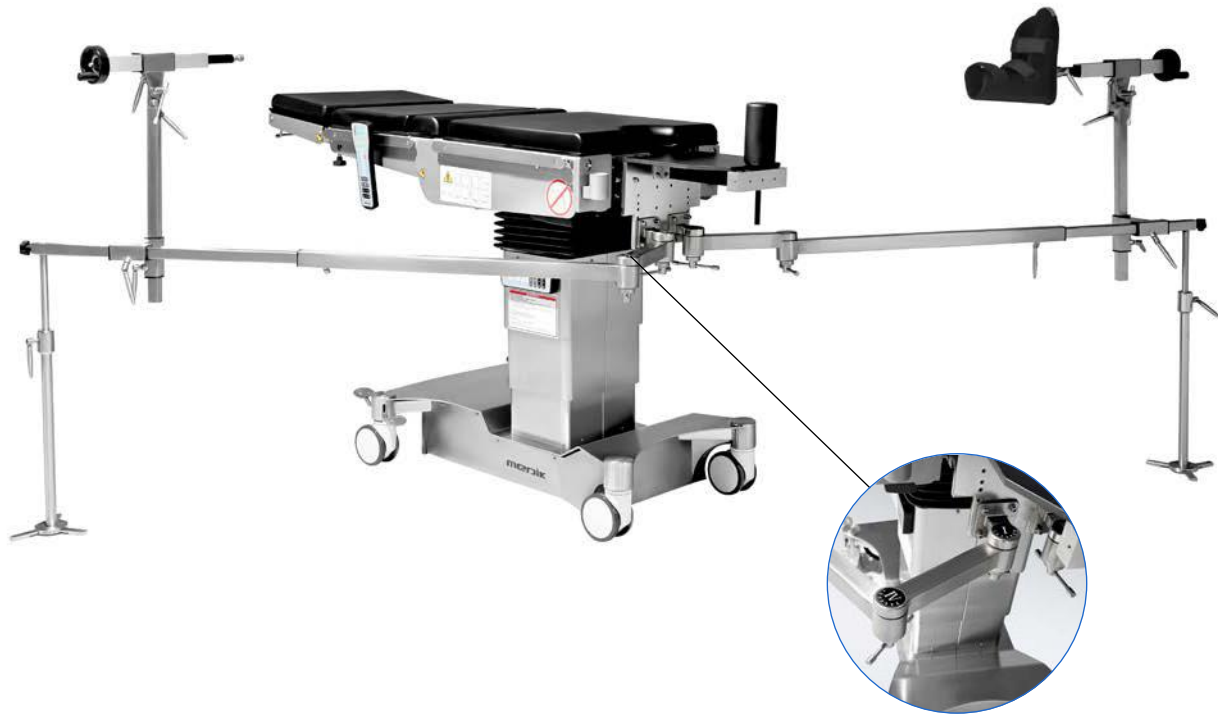

Seat board



Foot Adjustment Device

The foot section can be adjusted in multiple angles and distance, as well as the height, allowing doctors to easily fix the patient in the best surgical position.

PRODUCT DETAILS

Auxiliary trolley

Equipped with an auxiliary trolley for easy installation and use.



Traction crank

The traction crank can be used to traction and stretch the joints after fixing the patient's foot, With precise scale indication, making it easier for doctors to perform surgery.



Height adjustment device

The height of the traction frame can be adjusted according to the height of the operating table so that they remain balanced.

ABOUT US

Medik is a group focusing on surgical room solution and nursing care solution. We started to manufacture surgical table & lamp since 2011 with a professional team consisting of experienced designer and engineer from international top brand. With 10 years accumulation, Medik is becoming a new star to provide reliable product and valuable solution for surgical room.

Qualified with ISO9001, ISO13485, CE, Medik has successfully distributed the surgical products covering surgical table & lamp & pendant to international hospitals both domestic and overseas including many countries in Europe.



YOU TUBE



MEDIK

Pls scan the QR code for more info

www.mediksurgical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US